


CA[®] Pro / CA[®] Pro+

The powerful thermoforming material with flexible elastomer layer.

More endurance, less pressure: CA[®] Pro is setting new standards in aligner treatment. A flexible elastomer core in a hard-elastic double shell ensures a continuous force level at a reduced initial force and higher patient comfort.

CA[®] Pro can be used as single-splint aligner material as well as within the 3-splint system according to the CA[®] protocol.

- /// High elasticity and break resistance thanks to a flexible elastomer layer
- /// Continuous force transmission, minimal loss of power
- /// High wearing comfort thanks to gentle initial force
- /// Also available: CA[®] Pro+ with integrated protection foil

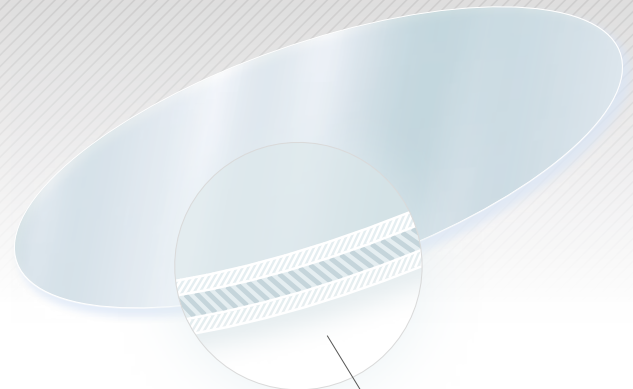

CA[®] Pro
Innovative double shell construction with elastomer core. Constant force level, less initial force.



Delivery program:

CA[®] Pro (without integrated protection foil)

/// CA [®] Pro	0.5 mm	Ø 125 mm	3640.1 10 pieces	3640.2 100 pieces
/// CA [®] Pro	0.625 mm	Ø 125 mm	3641.1 10 pieces	3641.2 100 pieces
/// CA [®] Pro	0.75 mm	Ø 125 mm	3642.1 10 pieces	3642.2 100 pieces

CE

CA[®] Pro (with integrated protection foil)

/// CA [®] Pro+	0.5 mm	Ø 125 mm	3644.1 10 pieces	3644.2 100 pieces
/// CA [®] Pro+	0.625 mm	Ø 125 mm	3645.1 10 pieces	3645.2 100 pieces
/// CA [®] Pro+	0.75 mm	Ø 125 mm	3646.1 10 pieces	3646.2 100 pieces

CE

CA[®] Pro+ is also available with a diameter of 120 mm.

IMPRIMO® LC Splint flex

The new 3D printer resin with thermoactive memory effect. Fits perfectly.


The new light-curing methacrylate-based 3D printer resin is indicated for the fabrication of high-precision, transparent occlusal splints and for bruxism and lower jaw protrusion splints. The thermoactive memory effect ensures perfect material fit.

When cured, the great flexibility provides high wearing comfort, low susceptibility to fracture and outstanding polishing properties.

Like almost all IMPRIMO® LC printer resins, IMPRIMO® LC Splint flex is also free of diluent monomers, enabling low-risk and safe handling for the user. Matched to DLP systems with a light source of 385 nm.

Features and benefits:

- /// High flexibility
- /// Thermo-active memory effect
- /// Excellent mechanical properties
- /// Parameters are available for Asiga MAX™, Asiga PRO 4K, RS wash, RS cure, IMPRIMO® Clean and IMPRIMO® Cure

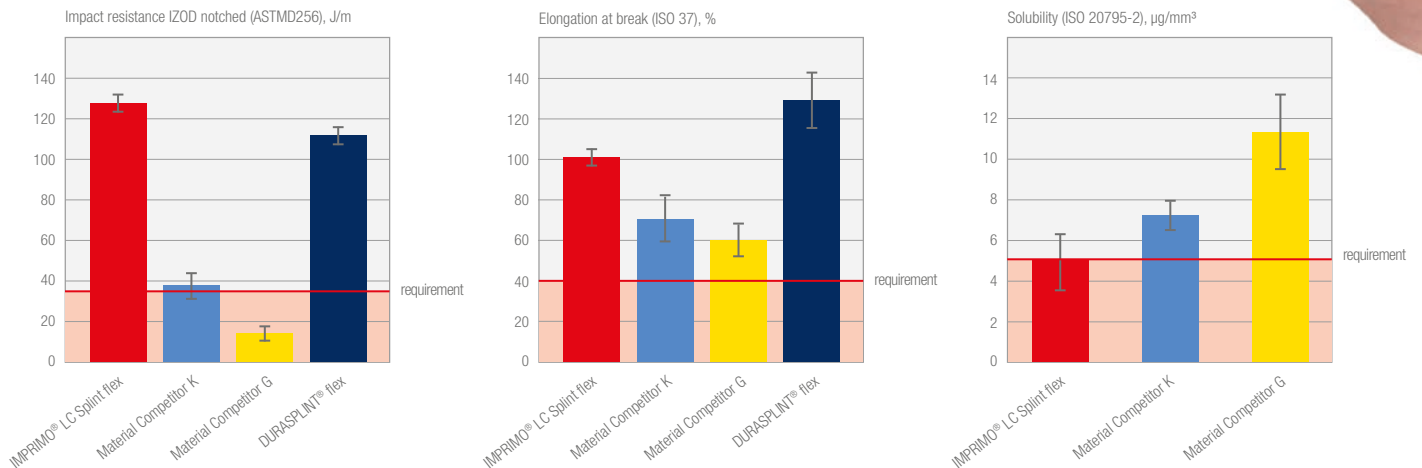

Delivery program:

/// IMPRIMO® LC Splint flex	500 g	6525.1
/// IMPRIMO® LC Splint flex	1000 g	6525.2

CE

Do you have questions? Visit our website www.scheu-dental.com/sales_export or get in touch with your local dealer.

IMPRIMO® Model Separator

The water-based, blue-transparent gel for insulating 3D-printed models. For an optimum insulating effect.

The IMPRIMO® Model Separator provides an optimum insulating effect between IMPRIMO® LC printer resins and STEADY RESIN autopolymers in the fabrication of removable orthodontic appliances. The insulating gel does not contain any volatile, organic solvents and thus ensures safe application. IMPRIMO® Model Separator is not suitable for 3D models generated by filament printer.

Features and benefits:

- //** Optimum insulating effect due to the formation of an even film on the model
- //** STL data set for 3D printing of the rack for pump bottles and brushes is available as download
- //** Available as 100 ml pump bottle with brush and as refill in the 500 ml refill bottle

Delivery program:

// IMPRIMO® Model Separator (pump bottle incl. brush)	100 ml	6551.1
// IMPRIMO® Model Separator REFILL (refill bottle)	500 ml	6551.2
// Brush Leonhardy	1 piece	6510.002



CA® TIP 5

Aligner pliers CA® TIP 5 for the correction of splint and aligner margins. Eliminate pressure points on the gingiva.

The aligner pliers CA® TIP 5 are specifically designed for adjusting aligner margins to eliminate any pressure points on the gingiva, making CA® TIP 5 ideal for the CA® PROFESSIONAL system, but also for any other systems incorporating gingiva-supported aligners.

Features and benefits:

- //** Ideal for finishing splint edges thanks to angled pliers head
- //** Sterilisable
- //** Made of stainless steel. 5-year warranty on material, processing, breakage and corrosion

Delivery program:

// CA® TIP 6	1 piece	5323.1
---------------------	---------	--------

CE




SCHEU

SCHEU GROUP

SMILE LITE FULL

CA® PROFESSIONAL FAMILY

CA
DIGITAL
SCHEUGROUP

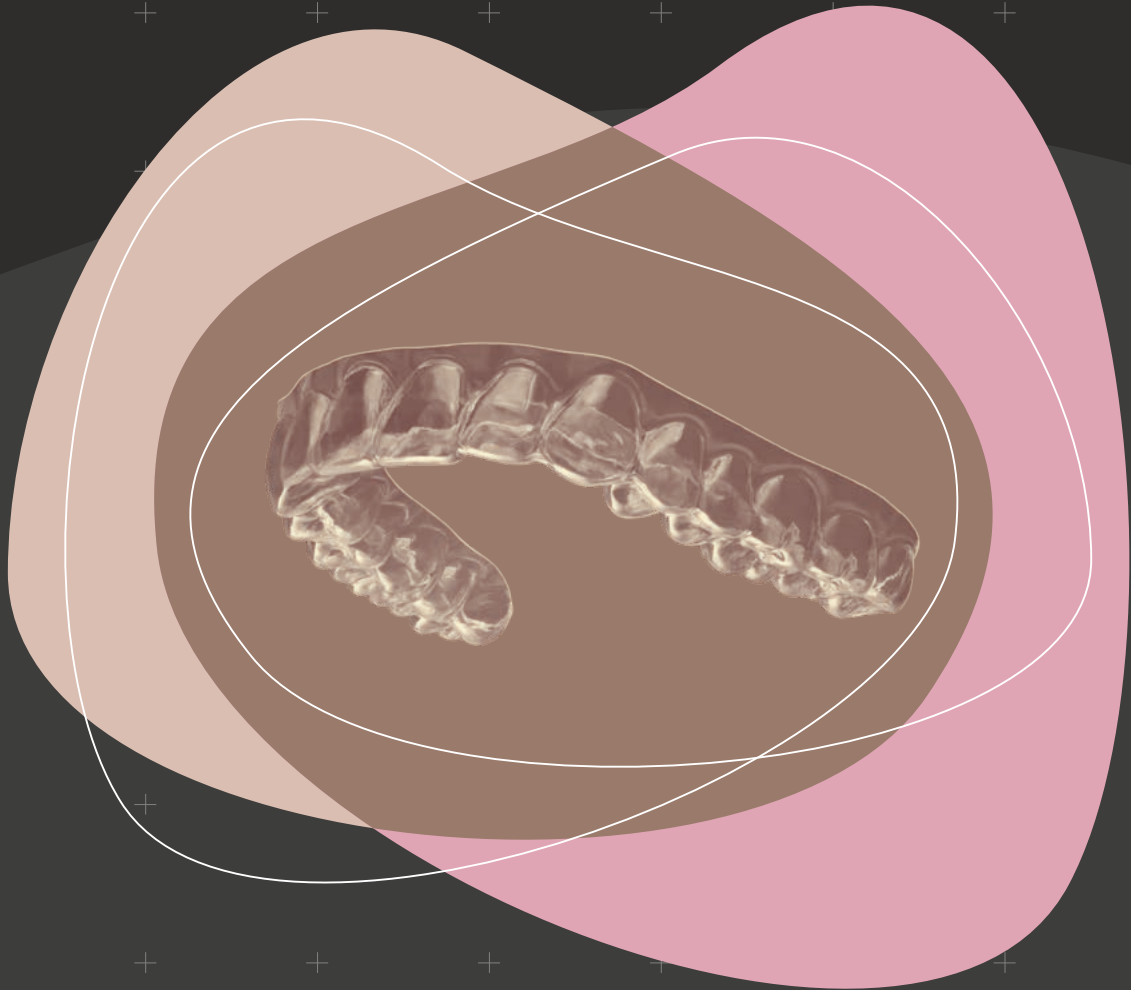

NOW FOR ALL INDICATIONS:
CA® PROFESSIONAL FAMILY

THE ALIGNER SYSTEM
WITH ALLROUND SUPPORT