

Safety Data Sheet

according to Regulation (EC) No 1907/2006

sili fluid Adhesive

Revision date: 13.02.2020

Product code: 10293_

Page 1 of 11

SECTION 1: Identification of the substance/mixture and of the company/undertaking

1.1. Product identifier

sili fluid Adhesive

1.2. Relevant identified uses of the substance or mixture and uses advised against

Use of the substance/mixture

Adhesive for use in dentistry.

1.3. Details of the supplier of the safety data sheet

Company name:	DETAX GmbH & Co. KG	
Street:	Carl-Zeiss-Strasse	
Place:	D-76275 Ettlingen	
Telephone:	+49 7243/510-0	Telefax: +49 7243/510-100
e-mail:	post@detax.de	
Internet:	www.detax.de	
Responsible Department:	Emergency number:	

+49 7243/510-0

This number is only obtainable during office hours (Monday - Thursday 8.00 a.m. - 5.00 p.m., Friday 8.00 a.m. - 4.00 p.m.)

1.4. Emergency telephone number:

+49 7243/510-0

This number is only obtainable during office hours (Monday - Thursday 8.00 a.m. - 5.00 p.m., Friday 8.00 - 4.00 p.m.)

SECTION 2: Hazards identification

2.1. Classification of the substance or mixture

Regulation (EC) No. 1272/2008

Hazard categories:

Flammable liquid: Flam. Liq. 2

Aspiration hazard: Asp. Tox. 1

Skin corrosion/irritation: Skin Irrit. 2

Serious eye damage/eye irritation: Eye Irrit. 2

Specific target organ toxicity - single exposure: STOT SE 3

Hazardous to the aquatic environment: Aquatic Chronic 2

Hazard Statements:

Highly flammable liquid and vapour.

May be fatal if swallowed and enters airways.

Causes skin irritation.

Causes serious eye irritation.

May cause drowsiness or dizziness.

Toxic to aquatic life with long lasting effects.

2.2. Label elements

Regulation (EC) No. 1272/2008

Hazard components for labelling

Naphtha (petroleum), hydrotreated light; Low boiling point hydrogen treated naphtha

ethyl acetate

methyl acetate

n-hexane

Signal word: Danger

Safety Data Sheet

according to Regulation (EC) No 1907/2006


sili fluid Adhesive

Revision date: 13.02.2020

Product code: 10293_

Page 2 of 11

Pictograms:

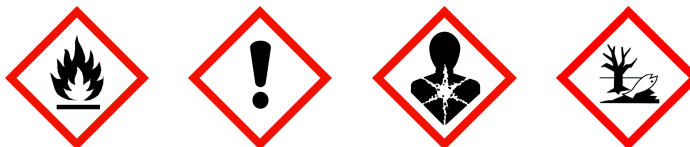

Hazard statements

H225	Highly flammable liquid and vapour.
H304	May be fatal if swallowed and enters airways.
H315	Causes skin irritation.
H319	Causes serious eye irritation.
H336	May cause drowsiness or dizziness.
H411	Toxic to aquatic life with long lasting effects.

Precautionary statements

P210	Keep away from heat, hot surfaces, sparks, open flames and other ignition sources. No smoking.
P273	Avoid release to the environment.
P301+P310	IF SWALLOWED: Immediately call a POISON CENTER/doctor.
P331	Do NOT induce vomiting.
P370+P378	In case of fire: Use Carbon dioxide (CO ₂), Foam, Extinguishing powder to extinguish.
P391	Collect spillage.
P403+P235	Store in a well-ventilated place. Keep cool.

Additional advice on labelling

According to Regulation (EC) 1272/2008, art.1 No. 5 (d) this product as a medical product must not be labelled!

2.3. Other hazards

No information available.

SECTION 3: Composition/information on ingredients

3.2. Mixtures

Chemical characterization

Alkyl silicone resin in organic solvents + pigment.

Safety Data Sheet

according to Regulation (EC) No 1907/2006

sili fluid Adhesive

Revision date: 13.02.2020

Product code: 10293_

Page 3 of 11

Hazardous components

CAS No	Chemical name			Quantity
	EC No	Index No	REACH No	
	GHS Classification			
64742-49-0	Naphtha (petroleum), hydrotreated light; Low boiling point hydrogen treated naphtha			25 - < 30 %
	265-151-9	649-328-00-1		
	Flam. Liq. 2, Skin Irrit. 2, STOT SE 3, Asp. Tox. 1, Aquatic Chronic 2; H225 H315 H336 H304 H411			
141-78-6	ethyl acetate			25 - < 30 %
	205-500-4	607-022-00-5	02-2119752482-38	
	Flam. Liq. 2, Eye Irrit. 2, STOT SE 3; H225 H319 H336 EUH066			
79-20-9	methyl acetate			5 - < 10 %
	201-185-2	607-021-00-X	02-2119752526-34	
	Flam. Liq. 2, Eye Irrit. 2, STOT SE 3; H225 H319 H336 EUH066			
110-54-3	n-hexane			1 - < 5 %
	203-777-6	601-037-00-0		
	Flam. Liq. 2, Repr. 2, Skin Irrit. 2, STOT SE 3, STOT RE 2, Asp. Tox. 1, Aquatic Chronic 2; H225 H361f H315 H336 H373 H304 H411			

Full text of H and EUH statements: see section 16.

SECTION 4: First aid measures

4.1. Description of first aid measures

After inhalation

Provide fresh air. In case of accident or unwellness, seek medical advice immediately (show directions for use or safety data sheet if possible).

After contact with skin

After contact with skin, wash immediately with plenty of water and soap. Take off immediately all contaminated clothing and wash it before reuse. Medical treatment necessary.

After contact with eyes

After contact with the eyes, rinse with water with the eyelids open for a sufficient length of time, then consult an ophthalmologist immediately.

After ingestion

Observe risk of aspiration if vomiting occurs. Rinse mouth immediately and drink plenty of water. Induce vomiting when the affected person is not unconscious. Medical treatment necessary. Seek immediately medical advice. Do not induce vomiting. In case of spontaneous vomiting take care of an unhindered flow out of the vomit (danger of suffocation).

4.2. Most important symptoms and effects, both acute and delayed

No information available.

4.3. Indication of any immediate medical attention and special treatment needed

Treat symptomatically.

SECTION 5: Firefighting measures

5.1. Extinguishing media

Suitable extinguishing media

Carbon dioxide (CO₂), Foam, Extinguishing powder.

according to Regulation (EC) No 1907/2006

sili fluid Adhesive

Revision date: 13.02.2020

Product code: 10293_

Page 4 of 11

Unsuitable extinguishing media

Water.

5.2. Special hazards arising from the substance or mixture

Highly flammable. Vapours can form explosive mixtures with air. Flammable. Vapours can form explosive mixtures with air.

5.3. Advice for firefighters

Wear a self-contained breathing apparatus and chemical protective clothing. Full protection suit.

Additional information

Use water spray jet to protect personnel and to cool endangered containers. Suppress gases/vapours/mists with water spray jet. Collect contaminated fire extinguishing water separately. Do not allow entering drains or surface water.

SECTION 6: Accidental release measures**6.1. Personal precautions, protective equipment and emergency procedures**

Remove all sources of ignition. Provide adequate ventilation. Do not breathe gas/fumes/vapour/spray. Avoid contact with skin, eyes and clothes. Use personal protection equipment. Remove all sources of ignition.

6.2. Environmental precautions

Do not allow uncontrolled discharge of product into the environment. Explosion risk in case of fire. Do not allow uncontrolled discharge of product into the environment..

6.3. Methods and material for containment and cleaning up

Absorb with liquid-binding material (e.g. sand, diatomaceous earth, acid- or universal binding agents). Treat the recovered material as prescribed in the section on waste disposal.

6.4. Reference to other sections

Safe handling: see section 7

Personal protection equipment: see section 8

Disposal: see section 13

SECTION 7: Handling and storage**7.1. Precautions for safe handling****Advice on safe handling**

If handled uncovered, arrangements with local exhaust ventilation have to be used. Do not breathe gas/fumes/vapour/spray.

Advice on protection against fire and explosion

Keep away from sources of ignition. - No smoking. Take precautionary measures against static discharges. Vapours can form explosive mixtures with air.

7.2. Conditions for safe storage, including any incompatibilities**Requirements for storage rooms and vessels**

Keep container tightly closed. Keep locked up. Store in a place accessible by authorized persons only. Provide adequate ventilation as well as local exhaust at critical locations. Keep in a cool, well-ventilated place. Keep away from heat, hot surfaces, sparks, open flames and other ignition sources. No smoking.

Hints on joint storage

Do not store together with: Pyrophoric or self-heating substances.

Do not store with combustible materials.

Further information on storage conditions

Keep only in the original container in a cool, dry and well-ventilated place, away from foodstuffs.

7.3. Specific end use(s)

Liquid for better adhesion of condensation curing dental impression materials on impression trays.

For use by trained specialist staff.

according to Regulation (EC) No 1907/2006

sili fluid Adhesive

Revision date: 13.02.2020

Product code: 10293_

Page 5 of 11

SECTION 8: Exposure controls/personal protection

8.1. Control parameters

Exposure limits (EH40)

CAS No	Substance	ppm	mg/m ³	fibres/ml	Category	Origin
141-78-6	Ethyl acetate	200	-		TWA (8 h)	WEL
		400	-		STEL (15 min)	WEL
79-20-9	Methyl acetate	200	616		TWA (8 h)	WEL
		250	770		STEL (15 min)	WEL
110-54-3	n-Hexane	20	72		TWA (8 h)	WEL

8.2. Exposure controls

Appropriate engineering controls

If handled uncovered, arrangements with local exhaust ventilation have to be used. Do not breathe gas/fumes/vapour/spray.

Protective and hygiene measures

Remove contaminated, saturated clothing immediately. Draw up and observe skin protection programme. Wash hands and face before breaks and after work and take a shower if necessary. When using do not eat, drink, smoke, sniff.

Eye/face protection

Suitable eye protection: goggles.

Hand protection

When handling with chemical substances, protective gloves must be worn with the CE-label including the four control digits. The quality of the protective gloves resistant to chemicals must be chosen as a function of the specific working place concentration and quantity of hazardous substances. For special purposes, it is recommended to check the resistance to chemicals of the protective gloves mentioned above together with the supplier of these gloves. Suitable are gloves of the following material: Butyl caoutchouc (butyl rubber), FKM (fluoro rubber)

Skin protection

Flame-retardant protective clothing. Wear anti-static footwear and clothing . .

Respiratory protection

In case of inadequate ventilation wear respiratory protection.

SECTION 9: Physical and chemical properties

9.1. Information on basic physical and chemical properties

Physical state: liquid:
 Colour: blue
 Odour: faintly like esters

Test method

pH-Value: not determined

Changes in the physical state

Melting point: not determined

Initial boiling point and boiling range: 65 °C DIN 51356

Flash point: -1 °C DIN 51755

Flammability

according to Regulation (EC) No 1907/2006

sili fluid Adhesive

Revision date: 13.02.2020

Product code: 10293_

Page 6 of 11

Solid: not applicable
 Gas: not applicable

Explosive properties

The product is not: Explosive.

Lower explosion limits: 1,5 vol. %
 Upper explosion limits: 12 vol. %
 Ignition temperature: 340 °C DIN 51794

Auto-ignition temperature

Solid: not applicable
 Gas: not applicable

Decomposition temperature: not determined

Oxidizing properties

Not oxidizing.

Vapour pressure: 135 hPa
 (at 20 °C)
 Vapour pressure: 490 hPa
 (at 50 °C)

Density (at 20 °C): 0,93 g/cm³ DIN 51757

Water solubility: partially miscible

Solubility in other solvents

not determined

Partition coefficient: not determined

Vapour density: not determined

Evaporation rate: not determined

9.2. Other information

Solid content: not determined

SECTION 10: Stability and reactivity

10.1. Reactivity

Highly flammable.

10.2. Chemical stability

The product is stable under storage at normal ambient temperatures.

10.3. Possibility of hazardous reactions

Reacts with : combustible materials.

10.4. Conditions to avoid

Keep away from sources of heat (e.g. hot surfaces), sparks and open flames. Vapours can form explosive mixtures with air. Higher temperatures advance the formulation of explosive vapour-air mixtures, therefore don't expose the product to increased temperatures.

10.5. Incompatible materials

No information available.

10.6. Hazardous decomposition products

The following applies for the silicone content of the product: At temperature of appr. 150°C/ 302 °F a small amount of formaldehyde can be released by oxidative degradation.

SECTION 11: Toxicological information

according to Regulation (EC) No 1907/2006

sili fluid Adhesive

Revision date: 13.02.2020

Product code: 10293_

Page 7 of 11

11.1. Information on toxicological effects

Acute toxicity

Based on available data, the classification criteria are not met.

CAS No	Chemical name					
	Exposure route	Dose	Species	Source	Method	
64742-49-0	Naphtha (petroleum), hydrotreated light; Low boiling point hydrogen treated naphtha					
	oral	LD50 mg/kg	2000	Rat		
	dermal	LD50 mg/kg	4000	Rabbit		
	inhalation (4 h) vapour	LC50	54 mg/l	Rat		
141-78-6	ethyl acetate					
	oral	LD50 mg/kg	5620	Rat	GESTIS	
	dermal	LD50 mg/kg	>18000	Rabbit	GESTIS	
79-20-9	methyl acetate					
	oral	LD50 mg/kg	>5000	Rat	GESTIS	
	dermal	LD50 mg/kg	>5000	Rabbit	GESTIS	
	inhalation (4 h) vapour	LC50	49-98	Rat		
110-54-3	n-hexane					
	oral	LD50 mg/kg	25000	Rat	GESTIS	
	inhalation (4 h) vapour	LC50	169 mg/l	Rat	GESTIS	

Irritation and corrosivity

Causes skin irritation.

Causes serious eye irritation.

Sensitising effects

Based on available data, the classification criteria are not met.

Carcinogenic/mutagenic/toxic effects for reproduction

Based on available data, the classification criteria are not met.

STOT-single exposure

May cause drowsiness or dizziness. (Naphtha (petroleum), hydrotreated light; Low boiling point hydrogen treated naphtha; ethyl acetate)

STOT-repeated exposure

Based on available data, the classification criteria are not met.

Aspiration hazard

May be fatal if swallowed and enters airways.

Additional information on tests

This mixture is classified as hazardous according to regulation (EC) No. 1272/2008 [CLP]. Special hazards arising from the substance or mixture!

SECTION 12: Ecological information

12.1. Toxicity

according to Regulation (EC) No 1907/2006

sili fluid Adhesive

Revision date: 13.02.2020

Product code: 10293_

Page 8 of 11

Toxic to aquatic life with long lasting effects.

CAS No	Chemical name					
	Aquatic toxicity	Dose	[h] [d]	Species	Source	Method
141-78-6	ethyl acetate					
	Acute fish toxicity	LC50 328 mg/l	96 h	Pimephales promelas (fathead minnow)	GESTIS	
	Acute algae toxicity	ErC50 mg/l 2500	96 h		GESTIS	
	Acute crustacea toxicity	EC50 679 mg/l	48 h	Daphnia magna (Big water flea)	GESTIS	
79-20-9	methyl acetate					
	Acute fish toxicity	LC50 320 mg/l	96 h	Pimephales promelas (fathead minnow)		
	Acute algae toxicity	ErC50 mg/l >120	72 h	Desmodesmus subspicatus.	OECD 201	
	Acute crustacea toxicity	EC50 mg/l 1027	48 h	Daphnia magna (Big water flea)	OECD 202	
110-54-3	n-hexane					
	Acute fish toxicity	LC50 mg/l 57,8	96 h	Pimephales promelas	GESTIS	

12.2. Persistence and degradability

The product has not been tested.

CAS No	Chemical name			
	Method	Value	d	Source
	Evaluation			
141-78-6	ethyl acetate			
	OECD 301D/ EEC 92/69/V, C.4-E	>70%	28	

12.3. Bioaccumulative potential

The product has not been tested.

Partition coefficient n-octanol/water

CAS No	Chemical name	Log Pow
110-54-3	n-hexane	3,9

12.4. Mobility in soil

The product has not been tested.

12.5. Results of PBT and vPvB assessment

Not identified as PBT/ vPvB substances

12.6. Other adverse effects

No information available.

Further information

Do not allow to enter into surface water or drains. Do not allow to enter into soil/subsoil.

SECTION 13: Disposal considerations

13.1. Waste treatment methods

Advice on disposal

Do not allow to enter into surface water or drains. Do not allow to enter into soil/subsoil. Dispose of waste according to applicable legislation.

according to Regulation (EC) No 1907/2006

sili fluid Adhesive

Revision date: 13.02.2020

Product code: 10293_

Page 9 of 11

Contaminated packaging

Non-contaminated packages may be recycled.

Handle contaminated packages in the same way as the substance itself.

SECTION 14: Transport information

Land transport (ADR/RID)

14.1. UN number:	UN 1866
14.2. UN proper shipping name:	Resin solution
14.3. Transport hazard class(es):	3
14.4. Packing group:	II
Hazard label:	3
Limited quantity:	5 L/ 30 kg
Transport category:	2
Hazard No:	33
Tunnel restriction code:	D/E

Other applicable information (land transport)

Flammable liquid

Marine transport (IMDG)

14.1. UN number:	UN 1866
14.2. UN proper shipping name:	Resin solution
14.3. Transport hazard class(es):	3
14.4. Packing group:	II
Hazard label:	3
Marine pollutant:	-
Special Provisions:	-
Limited quantity:	5 L/ 30 kg
EmS:	F-E, S-E

Other applicable information (marine transport)

Flash point: -1°C c.c.

Air transport (ICAO-TI/IATA-DGR)

14.1. UN number:	UN 1866
14.2. UN proper shipping name:	Resin solution
14.3. Transport hazard class(es):	3
14.4. Packing group:	II
Hazard label:	3
Limited quantity Passenger:	1 L/ 30 kg
Passenger LQ:	Y341
Excepted quantity:	E2
IATA-packing instructions - Passenger:	353
IATA-max. quantity - Passenger:	5 L
IATA-packing instructions - Cargo:	364
IATA-max. quantity - Cargo:	60 L

Other applicable information (air transport)

Flammable liquid

14.6. Special precautions for user

Warning: Combustible liquid.

according to Regulation (EC) No 1907/2006

sili fluid Adhesive

Revision date: 13.02.2020

Product code: 10293_

Page 10 of 11

14.7. Transport in bulk according to Annex II of Marpol and the IBC Code

not applicable

SECTION 15: Regulatory information**15.1. Safety, health and environmental regulations/legislation specific for the substance or mixture****EU regulatory information**

Information according to 2012/18/EU (SEVESO III): E2 Hazardous to the Aquatic Environment

Additional information: P5c

National regulatory information

Employment restrictions: Observe restrictions to employment for juvenils according to the 'juvenile work protection guideline' (94/33/EC).

Water contaminating class (D): 2 - clearly water contaminating

15.2. Chemical safety assessment

Chemical safety assessments for substances in this mixture were not carried out.

SECTION 16: Other information**Abbreviations and acronyms**

ADR: Accord européen sur le transport des marchandises dangereuses par Route (European Agreement concerning the International Carriage of Dangerous Goods by Road)

IMDG: International Maritime Code for Dangerous Goods

IATA: International Air Transport Association

GHS: Globally Harmonized System of Classification and Labelling of Chemicals

EINECS: European Inventory of Existing Commercial Chemical Substances

ELINCS: European List of Notified Chemical Substances

CAS: Chemical Abstracts Service

LC50: Lethal concentration, 50%

LD50: Lethal dose, 50%

CLP: Classification, labelling and Packaging

REACH: Registration, Evaluation and Authorization of Chemicals

GHS: Globally Harmonised System of Classification, Labelling and Packaging of Chemicals

UN: United Nations

DNEL: Derived No Effect Level

DMEL: Derived Minimal Effect Level

PNEC: Predicted No Effect Concentration

ATE: Acute toxicity estimate

LL50: Lethal loading, 50%

EL50: Effect loading, 50%

EC50: Effective Concentration 50%

ErC50: Effective Concentration 50%, growth rate

NOEC: No Observed Effect Concentration

BCF: Bio-concentration factor

PBT: persistent, bioaccumulative, toxic

vPvB: very persistent, very bioaccumulative

RID: Regulations concerning the international carriage of dangerous goods by rail

EmS: Emergency Schedules

MFAG: Medical First Aid Guide

ICAO: International Civil Aviation Organization

MARPOL: International Convention for the Prevention of Marine Pollution from Ships

IBC: Intermediate Bulk Container

Safety Data Sheet

according to Regulation (EC) No 1907/2006

sili fluid Adhesive

Revision date: 13.02.2020

Product code: 10293_

Page 11 of 11

VOC: Volatile Organic Compounds
 SVHC: Substance of Very High Concern
 @1602.B016012

Classification for mixtures and used evaluation method according to Regulation (EC) No. 1272/2008 [CLP]

Classification	Classification procedure
Flam. Liq. 2; H225	On basis of test data
Asp. Tox. 1; H304	Calculation method
Skin Irrit. 2; H315	Calculation method
Eye Irrit. 2; H319	Calculation method
STOT SE 3; H336	Calculation method
Aquatic Chronic 2; H411	Calculation method

Relevant H and EUH statements (number and full text)

H225	Highly flammable liquid and vapour.
H304	May be fatal if swallowed and enters airways.
H315	Causes skin irritation.
H319	Causes serious eye irritation.
H336	May cause drowsiness or dizziness.
H361f	Suspected of damaging fertility.
H373	May cause damage to organs through prolonged or repeated exposure.
H411	Toxic to aquatic life with long lasting effects.
EUH066	Repeated exposure may cause skin dryness or cracking.

Further Information

The information is based on present level of our knowledge. It does not, however, give assurances of product properties and establishes no contract legal rights. The receiver of our product is singularly responsible for adhering to existing laws and regulations.

(The data for the hazardous ingredients were taken respectively from the last version of the sub-contractor's safety data sheet.)